EC1260BL

Combination Warmer

The EC1260BL is designed to safely warm and store blankets and either irrigation fluids or injection fluids.

The Dual-chambered warming cabinet is constructed with stainless steel exterior casing and door with handle and hinges designed to withstand heavy usage. The cabinet is warmed using a low-heat-density electrothermal cable array. The electrothermal cable is positioned in the floor and two sides of the warming cabinet, providing even heating of the interior chamber.

WarmRight: Provides the right temperature for the right product.

Durable: The affordable, rugged stainless steel cabinets are made to last and work just as hard as you do.

Glass Door: Energy-efficient coated glass door reflects heat back into the cavity.

Service: Our worldwide support team is dedicated to providing quality equipment. We walk with you every step of the way to maximize your return on investment.

Made in the USA with a commitment to quality

An ISO 13485:2016 certified company

Standard Features

- Windowed doors
- Epoxy-coated interior surfaces
- Programmable control
- Timer feature
- Adjustable interior lighting
- Temperature management access point
- Audible over-temp alarm
- Four (4) leveling feet



28

18 (1) liter bottles or 28 (1) liter bags

30

25-30 Blankets





Adjustable temperature range of

37°-40°C (98°-104°F) injection mode 37°-66°C (98°-150°F) irrigation mode 37°-93°C (98°-200°F) blanket cavity

Accurate temperature

IRR +0/-1.67°C (-3°F)

INJ +0/-1.10°C (-2°F)





Configurations (select one)

Door Assembly

■ Widowed door (standard)
■ Solid door

Door Swing

☐ Right hinged (standard) ☐ Left hinged

Door Lock

☐ No lock (standard) ☐ Manual lock

Electrical

□ 120V 1PH □ 230V 1PH

Accessories (select all that apply)

Seismic Mounting Brackets

Special installation requirements apply


■ Contact factory


Built in Trim Kit


☐ Trim Kit)



Specification

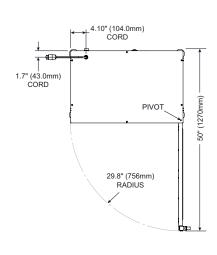

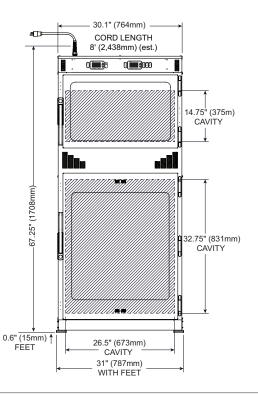

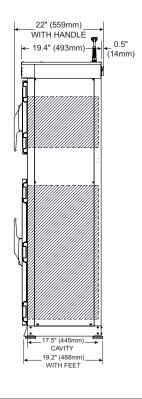



CLEARANCE

Top: 1"

Left: 1"

Right: 1"

Back: 3"



Model Exterior (H x W x D) Interior (H x W x D)

Net Weight

EC1260BL

67.25"(1708mm) x 31"(787mm) x 22"(559mm)

Ship Dimensions (L x W x H)*

47.5"(1207mm) x 26.5"(673mm) x 17.25"(438mm)

Dl.... C - ... fi -

320 lbs (145 kg)

DIMENSIONS

Ship Weight*

38" x 35"x 80" 420 lbs (191 kg) EC1260BL

*Domestic ground shipping information. Contact factory for export weight and dimensions.



	V	Pn	HZ	А	KVV	Plug Configuration
EC1260BL	120	1	50/60	14.2	1.7	NEMA 5-20P; 20A 125V Plug; Hospital Grade
	230	1	50/60	6.9	1.58	BS 1363 Plug; CEE 7/7 220-230V Plug



*Other international plugs are available, contact the manufacturer for more information.

**Electrical connections must meet all applicable federal, state and local codes.

IPX-0