

EC1540BL

Combination Warmer



The EC1540BL is designed to safely warm and store blankets and either irrigation fluids or injection fluids.

The dual-chambered warming cabinet is constructed with stainless steel exterior casing and door with handle and hinges designed to withstand heavy usage. The cabinet is warmed using a low-heat-density electrothermal cable array. The electrothermal cable is positioned in the floor and two sides of the warming cabinet, providing even heating of the interior chamber.

WarmRight: Provides the right temperature for the right product.

Durable: The affordable, rugged stainless steel cabinets are made to last and work just as hard as you do.

Glass Door: Energy-efficient coated glass door reflects heat back into the cavity.

Mobile: Rugged casters are designed to easily transport warmed fluids and blankets throughout the facility.

Service: Our worldwide support team is dedicated to providing quality equipment. We walk with you every step of the way to maximize your return on investment.

Made in the USA with a commitment to quality

An ISO 13485:2016 certified company

Standard Features

- Windowed door
- Organizational fluid baskets
- Epoxy-coated interior surfaces
- Full perimeter bumper
- Programmable control
- Timer feature
- Temperature management access point
- Audible over-temp alarm
- Four casters, 2 with brake



72 54 (1) liter bottles or 72 (1) liter bags

25 20-25 Blankets



Adjustable temperature range of

37°-40°C (98°-104°F) injection mode

37°-66°C (98°-150°F) irrigation mode

37°-93°C (98°-200°F) blanket cavity

Accurate temperature

IRR +0/-1.67°C (-3°F)

INJ +0/-1.10°C (-2°F)



Configurations (select one)

Door Assembly

- Window door (standard) Solid door

Door Swing

- Right hinged (standard) Left hinged

Door Lock

- No lock (standard) Manual lock

Electrical

- 120V 1PH 230V 1PH

Accessories (select all that apply)

Legs

- Leg kit, 6" (152mm) (5011149)

Seismic Mounting Brackets

Special installation requirements apply

- Contact factory


Built in Trim Kit

- Trim Kit

EC1540BL Combination Warming Cabinet



Specification

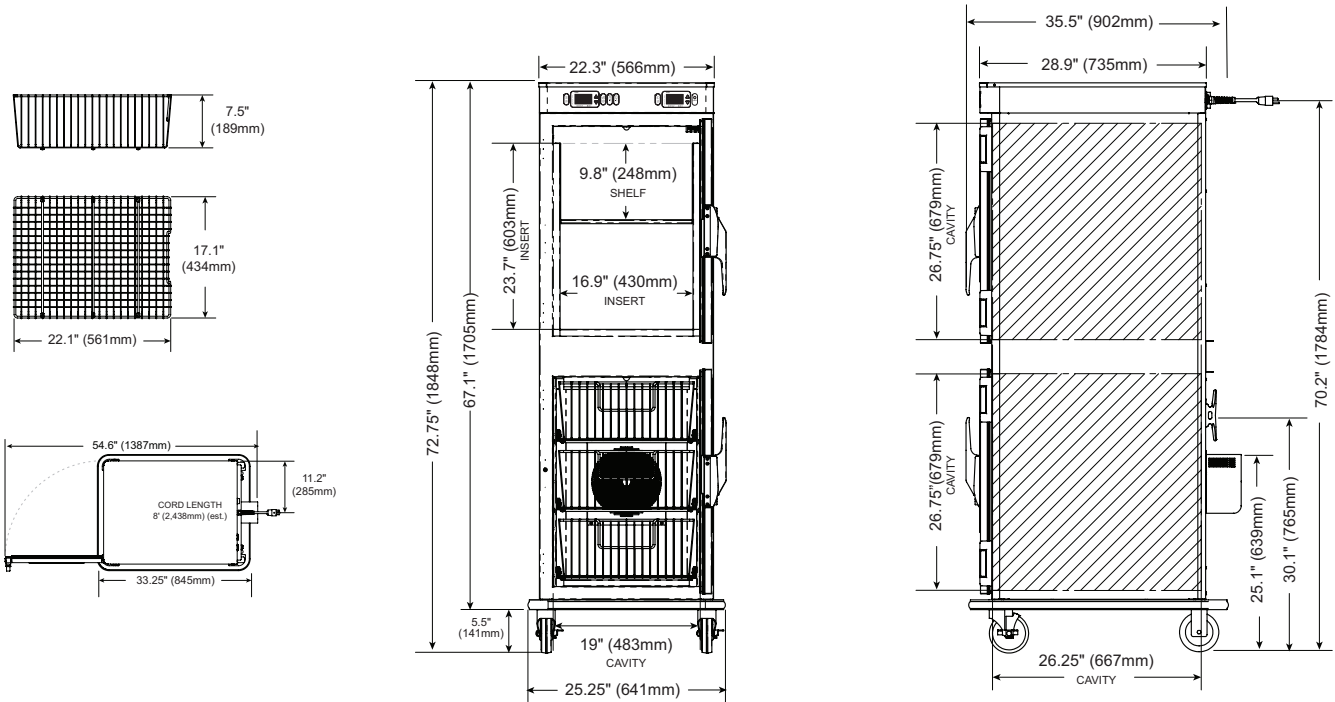

Top: 1"
Left: 1"
Right: 1"
Back: 3"



Model	Exterior (H x W x D)	Interior (H x W x D)	Net Weight
EC1540BL	72.75"(1848mm) x 24"(610mm) x 32.5"(826mm)	53.5"(1359mm) x 19"(483mm) x 26.25"(667mm)	374 lbs (170 kg)

	Ship Dimensions (L x W x H)*	Ship Weight*
EC1540BL	38" x 35" x 80"	494 lbs (224 kg)

*Domestic ground shipping information. Contact factory for export weight and dimensions.



	V	Ph	Hz	A	kW	Plug Configuration
EC1540BL	120	1	50/60	16	1.92	NEMA 5-20P; 20A 125V Plug; Hospital Grade
	230	1	50/60	7.4	1.7	BS 1363 Plug; CEE 7/7 220-230V Plug

*Other international plugs are available, contact the manufacturer for more information.

**Electrical connections must meet all applicable federal, state and local codes.



IPX-0

CONTACT US

W164 N9221 Water Street | Menomonee Falls, Wisconsin 53051 | U.S.A.
Phone: 262.251.8356 | 800.862.9276 | Fax: 262.251.7067 | enthermics.com